
Carnival Tube #158

USER'S GUIDE



For Technical Support:

Call our Technical Service Department
Monday through Friday, 9 a.m. to 5 p.m. (EST)
1-800-832-8697
customer_support@enablingdevices.com



50 Broadway
Hawthorne, NY 10532
Tel. 914.747.3070 / Fax 914.747.3480
Toll Free 800.832.8697
www.enablingdevices.com

Rainbows of colors!

A simple swipe at the mounted tube creates lights, music and vibration. It's a magical combination of sound and color that will captivate the young and old. Watch and listen as a rainbow of multi-colored beads and lights cascade through the tube. Great for encouraging reaching, swiping and grasping. Can also be used as a switch. It has a toggle switch which allows you to turn on or off any of the rewards. Size: 14"L x 6¾"W x 4½"H. Requires 4 AA Batteries. Weight: 1¾ lbs.

Operation:

1. The Carnival Tube requires 4 AA batteries (Not included). **Please Note:** Alkaline batteries are required (e.g. Duracell or Energizer brand) Do not use rechargeable or heavy duty batteries as they will not work with this device. Never mix old and new batteries together. Carefully turn the unit over, and then remove the Philips head screws from the battery covers. Install new batteries, being careful to observe proper battery polarity. Replace covers and screws.
2. To turn the desired stimuli ON or OFF (1-vibration; 2-music; 3-lights), push the small white rocker switches located on the bottom of the unit, to the ON or (open position). Open position means OFF. Please see Fig 1.

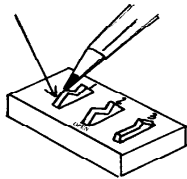


Fig .1

3. Rotate or spin the Carnival Tube and you will be rewarded with lights, music and vibration (if selected). Stopping rotation or spinning of the Carnival Tube will cause the Lights, Music and Vibration to stop after a brief delay of about 5 seconds.
4. To use the Carnival Tube as a switch plug one end of the supplied 1/8" to 1/8" cord into the Textured Roller's input jack

and the other end into your adapted toy/device. Rotate or spin the Carnival Tube and your toy/device will activate as long as rotating or spinning of the Roller is sustained.

Troubleshooting:

Problem: The Carnival Tube has no light or Vibration.

Action #1: Check both sets of batteries for strength and proper battery polarity. Replace if weak or dead.

Problem: The Carnival Tube has no music.

Action: Check both sets of batteries for strength and proper battery polarity. Replace if weak or dead.

Problem: The Carnival Tube does not work at all.

Action: Check all Batteries for strength and proper battery polarity. Replace if weak or dead.

Problem: Carnival Tube does not activate a toy or device.

Action # 1: Make sure all connections are tight.

Action # 2: Check batteries in toy/device, replace if weak or dead.

Action # 3: If adapters are being used they must be Mono, stereo adapters will not work.

Care of Unit:

The Carnival Tube can be wiped clean with mild household multi-purpose cleaner and disinfectant. We recommend Simple Green, which is a non toxic biodegradable all-purpose cleaner.

Do not submerge the unit in any fluids, as it will damage the contents and the electrical components.

Do not use abrasive cleaners, as they will scratch the surface of the unit.