

Remote Control Bubble Tube



Catalog No. 2229

enabling
devices

50 Broadway
Hawthorne, NY 10532
Tel. 914.747.3070 / Fax 914.747.3480
Toll Free 800.832.8697
www.enablingdevices.com

Changing colors and moving bubbles create a striking focal point for any room!

Remote control allows you to easily change colors and light patterns up to 15 feet away. Bubble speed is also adjustable. The Bubble Tube can provide either stimulation or promote relaxation depending upon the individual and the desired outcome. It also improves visual development, visual perception, color recognition and communication. *Size: 17"L x 17"W x 29"H. Uses AC Wall Transformer. Weight: 10 lbs.*

OPERATION:

Setup:

You will first need to place your Bubble Tube in the location where it will be used; these factors must be present for proper operation. One the Bubble Tube must be placed on a flat secure level surface, e.g. floor or very sturdy table, do not place the Bubble Tube on any type of folding table or unstable, un-level surface. The Bubble Tube weighs 30 Lbs. when filled with approximately 2.25 gallons of (Distilled water). Second the Bubble Tube requires a 120 Volt wall outlet AC power source. **Please Note:** If placing the Bubble Tube on the floor the power cord and any extension cords if used should be taped down to prevent any of its users from tripping on the cords.

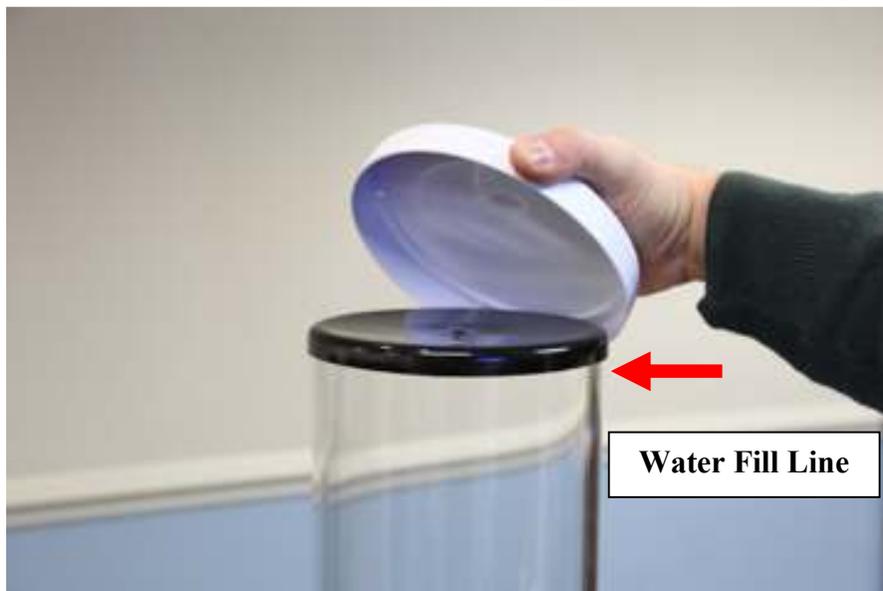


Photo No.1



Photo No.2

1. Now that you have placed your Bubble Tube in its location you will need to fill it with **(Distilled Water) Please Do Not Use Tap or Spring Water, doing so will cause algae to form inside the Bubble Tube and on its bubble aerator, causing a diminish in bubbles.** There are two caps that must first be removed before filling the Bubble Tube (See photos 1 & 2.). Lift off the first hard plastic cap; next remove the soft vinyl rubber sealing cap. You will need approximately 2.25 gallons of Distilled water to fill the Bubble Tube. Fill the tube up to 1 inch from the bottom of back cap. Do not fill past this, doing so will cause the water to bubble out of the tube. Once you have finished filling the Bubble Tube replace the two caps with the soft vinyl cap being put on first next the hard cap.



Photo No.3

To turn the Bubble Tube on flip the ON/OFF toggle switch on the base to ON. This will turn on the bubbles. The included IR remote control will let you tune the lights ON, OFF, change their color, brightness and how they light, Flash, Strobe, Fade and Smooth. The remote is powered by a CR2025 Watch battery which is included. Please remove the clear plastic battery tab before use.

Please Note: The remote is IR, meaning it is infrared just like a TV remote so it works in line of site. On the base of the Bubble Tube you will see a small depression with a small black IR receiver. You must aim your remote at this to control the Bubble Tube.

Bubble Adjustment:

The Bubble Adjustment knob allows you to control the intensity of the bubbles. Turning it clockwise increases the bubble flow, counter-clockwise decreases the bubble flow.

Draining the Bubble Tube:

If the Bubble Tube needs to be moved or you want to change the water, you will need to drain the water out. We have included a hand siphon pump to make this job easier. First you will need to remove the two top covers on the Bubble Tube and have a bucket or container capable of holding at least 2.5 gallons of water to drain the Bubble Tube into. The bucket or container must be placed below the Bubble Tube for the siphon to work. Next place the end of the siphon tube which is on the opposite side of the arrow indicator on the squeezing bulb (shown in photo below); inside the Bubble Tube. Make sure the tube goes all the way to the bottom of the Bubble Tube. Place the other end of the siphon tube into your bucket or sink and start the siphoning by squeezing the bulb until the water starts to flow out of the siphon tube into the bucket. Once this happens you can stop squeezing on the bulb and let the Bubble Tube drain out.



Troubleshooting:

If the unit is not functioning as it is supposed to:

- Check the AC power cord is plugged in the outlet all the way. Try unplugging the cord and then re-plug it back in.
- Make sure that the AC wall outlet is working.
- Check all your settings.
- If you have tried ALL of the above contact our Tech Support Department at ext. 320.



For Technical Support:

Call our Technical Service Department
Monday through Friday, 9 a.m. to 5 p.m. (EST)
1-800-832-8697
customer_support@enablingdevices.com

Rev 4/7/17