

Soft and sparkly!

Many therapists and teachers have been asking us to make a squeezable “Gooshy Switch”. We found a durable, soft, sparkly pillow and converted it into two switches: Gooshy Switch and a Gooshy Switch with Lights, Vibration, and Music. The Lights, Vibration and Music Switch has a special control which allows the user to select any or all of the rewards. It can be used alone as a sensory device. Size: 10"L x 4½"W x 2"H. Requires 2 AA Batteries. Weight: 1¼ lb.

Operation:

1. The Gooshey Switch requires two AA batteries. The battery compartment is located under the base of the unit. Carefully turn the switch over, and then remove the battery cover with a small Philips screwdriver. Install new batteries, being careful to observe proper battery polarity. Alkaline batteries are recommended because they last longer.
2. Replace the battery case cover and secure it with the small screw. Do not over-tighten.
3. To turn on or off the desired stimuli (1-vibration; 2-music; 3-lights), push the small white rocker switches to the on or (open position). Open position means off. Please see Fig 1.

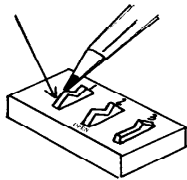


Figure 1.

4. Press anywhere on the Gooshey surface and you will be rewarded with the selected stimuli vibration, lights & music.

5. To operate as a switch connect an external toy or device via the 1/8-inch cord jack on the unit. If you need to use a ¼" adaptor, it must be a mono adaptor, not stereo. Press anywhere on the Gooshey surface to activate your toy or device. Once you release pressure from the switch your toy or device will turn off.

Troubleshooting:

Problem: The Gooshey Switch does not operate properly.

Action # 1: Make sure the batteries are in the battery compartment properly, and making good contact.

Action # 2: Make sure you have the desired stimuli selected as explained in operation No.3.

Action # 3: If the unit still does not work properly, replace the batteries.

Problem: The Gooshey Switch does not activate a connected toy/device.

Action #1: Make sure all connections are plugged in all the way. There should be no gaps.

Action #2: Check batteries in toy/device, replace if weak or dead.

Care of Unit:

The Gooshey Switch can be wiped clean with any household multi-purpose cleaner and disinfectant.

Do not submerge the unit, as it will damage the contents of the top and the electrical components.

Gooshey Switch #636

(With Vibration, Lights & Music)

USER'S GUIDE



For Technical Support:

Call our Technical Service Department
Monday through Friday, 9 a.m. to 5 p.m. (EST)
1-800-832-8697
customer_support@enablingdevices.com



50 Broadway
Hawthorne, NY 10532
Tel. 914.747.3070 / Fax 914.747.3480
Toll Free 800.832.8697
www.enablingdevices.com