
Multisensory Funny Face Switch

Works as a switch too!

Designed to provide auditory, tactile and visual stimulation. There are lights, music, and vibration when you push your tactile switch plate. Our Funny Face has a control on the bottom that allows you to turn any of the rewards on or off. Great cause and effect device. Can also be used as a capability switch to activate other devices. Size: 11"D. Requires 2 AA Batteries. Weight: 1½ lbs.

Operation:

1. The Funny Face Switch requires two AA batteries. The battery compartment is located under the base of the unit. Carefully turn the switch over, and then remove the battery cover with a small Philips screwdriver. Install new batteries, being careful to observe proper battery polarity. Alkaline batteries are recommended because they last longer.
2. Replace the battery case cover and secure it with the small screw. Do not over-tighten.
3. To turn on or off the desired stimuli (1-vibration; 2-music; 3-lights), push the small white rocker switches to the on or (open position). Open position means off. Please see Fig 1.

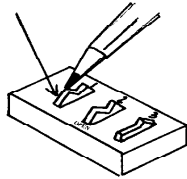


Figure 1.

4. Press anywhere on the Yellow textured plate and you will be rewarded with the selected stimuli vibration, lights & music.

5. To operate as a switch connect an external toy or device via the 1/8-inch cord jack on the unit. If you need to use a ¼" adaptor, it must be a mono adaptor, not stereo. Press anywhere on the textured yellow plate to activate your toy or device. Once you release pressure from the switch your toy or device will turn off.

Troubleshooting:

Problem: The Funny Face Switch does not operate properly.

Action # 1: Make sure the batteries are in the battery compartment properly, and making good contact.

Action # 2: Make sure you have the desired stimuli selected as explained in operation No.3.

Action # 3: If the unit still does not work properly, replace the batteries.

Problem: The Funny Face Switch does not activate a connected toy/device.

Action #1: Make sure all connections are tight.

Action #2: Check batteries in toy/device, replace if weak or dead.

Care of Unit:

The Funny Face Switch can be wiped clean with any household multi-purpose cleaner and disinfectant.

Do not submerge the unit, as it will damage the contents of the top and the electrical components.

Multisensory Funny Face Switch #743

(With Vibration, Lights & Music)

USER'S GUIDE



For Technical Support:

Call our Technical Service Department
Monday through Friday, 9 a.m. to 5 p.m. (EST)
1-800-832-8697

customer_support@enablingdevices.com



50 Broadway
Hawthorne, NY 10532
Tel. 914.747.3070 / Fax 914.747.3480
Toll Free 800.832.8697
www.enablingdevices.com