

Gentle touch, merry reward!

When the rim on the dome is gently pressed, silvery glitter swirls around while music, lights and vibration reward the user and provide them with auditory, visual and tactile stimulation. The sensory stimuli can be switched on or off individually. This dome also acts as a switch.

Size: 8"W x 4½"H. Requires 2 AA Batteries. Weight: 1¼ lbs.

OPERATING INSTRUCTIONS:

1. Carefully turn the unit over to reveal the battery compartment. Using a small Phillips-head screwdriver, remove the cover and install four "AA" batteries, being careful to observe proper polarity.
2. To turn on or off the desired stimuli (1 – vibration; 2 – music; 3 - lights), push the small white rocker switches to the on or (open position). Open position means off. Please see **Figure 1**.

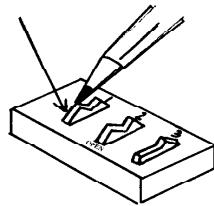


Figure 1

3. To use this unit as a switch to control other devices, plug the supplied cord into the jack on the back of the unit and connect the other end of the cord into the device you wish to operate. If you leave the stimuli on, they will operate simultaneously with the external device.

IMPORTANT NOTES:

- If the lights, vibration or motor stop working, please replace the batteries (4 AA cells).
- Make sure there are no obstacles between the clear dome and yellow base.
- When connecting a toy/device, be sure there are no gaps between any of the connections. All should be flush.
- Try your toy/device with a different switch to rule out the Saucer Dome Switch as the source of the problem.

If you have any problems or questions about this unit or any of our products, please call our **Technical Assistance Department** at:

(1-800-832-8697)

**Saucer Dome Switch w/Lights,
Vibration & Music #78**

User's Guide



For Technical Support:

Call our Technical Service Department
Monday through Friday, 9 a.m. to 5 p.m. (EST)
1-800-832-8697

customer_support@enablingdevices.com



50 Broadway
Hawthorne, NY 10532
Tel. 914.747.3070 / Fax 914.747.3480
Toll Free 800.832.8697
www.enablingdevices.com